

**Labels Specifications External Telescope Simulator 2.8X-FF Light**  
**( For EU -English, German, French, and Spanish)**


**1. RM00730-02 CE mark + European Authorized Representative label Rev.1**

- 1.1 Size of Label: 150\*31 mm ± 0.1 mm
- 1.2 Label Background Color: Black
- 1.3 Materials: Black Mat Vinyl; glue 851
- 1.4 European Rep. Printing color - White; Font type - Arial; Font size - as per illustration
- 1.5 CE - Font type "CE"; Font size see drawing (referenced dimensions)
- 1.6 Raw material of sticker: RM00728-02
- 1.7 UDI position: as per illustration



**2. RM00843-02 ETS 2.8X-FF Light Sticker on the Transportation Box (EU) Rev.2**

- 2.1 Size of Label: 150\*31 mm ± 0.1 mm
- 2.2 Material: Black Mat Vinyl; glue 851
- 2.3 Manufacturer printing color - White; Font - Arial; Font size -as pre illustration
- 2.4 Raw material of sticker: RM00728-02
- 2.5 UDI position: as per illustration

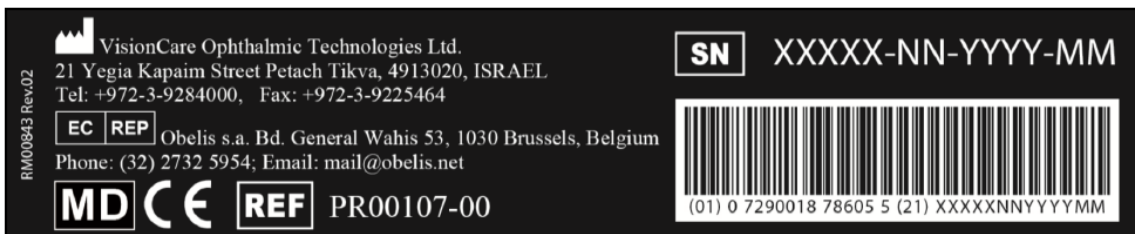

**Title: Labels Specifications for ETS 2.8X-FF Light (EU)**

Reference code: 100120006

Revision.0


**3. RM00874-02ETS 2.8X-FF Light Case Identification Label EU Rev.3**

- 3.1 Size of Label: 100\*120 mm ± 0.1 mm
- 3.2 Material: Polycarbonate (Lexan); Dual sided glue 3M 467
- 3.3 Thickness: 0.250mm
- 3.4 Logo: according to Specification No. 100050137-Specification of VisionCare Logo.
- 3.5 Label inscriptions: Color: print in color, Font Type: Ariel (6th row in Bold). Font size: as per illustration, Background: White, Label position: centered.



**4. RM01095-02 ETS Identification Label on telescope barrel (EU) Rev. 0**

- 4.1 Size of Label: 84\*14 mm ± 0.1 mm
- 4.2 Font type: Arial narrow
- 4.3 CE - Font type "CE"; Font size as per illustration
- 4.4 Raw material of sticker: RM01094-02
- 4.5 UDI position: as per illustration

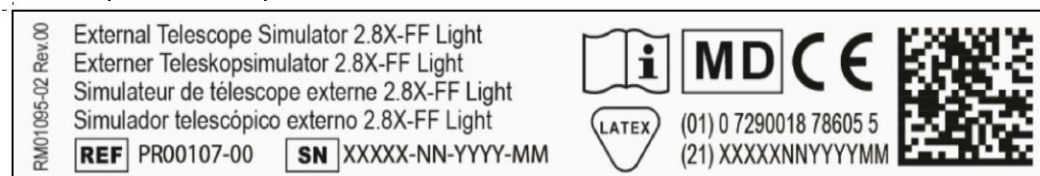

**Title: Labels Specifications for ETS 2.8X-FF Light (EU)**


Reference code: 100120006

Revision.0

**5. RM00274-02 FRAGILE LABEL Rev .0**



ITEM NO.	PartNo	DESCRIPTION	CatNo	QTY.
13	RD006-0106	ETS Identification Label on telescope barrel (EU)	RM01095-02	1

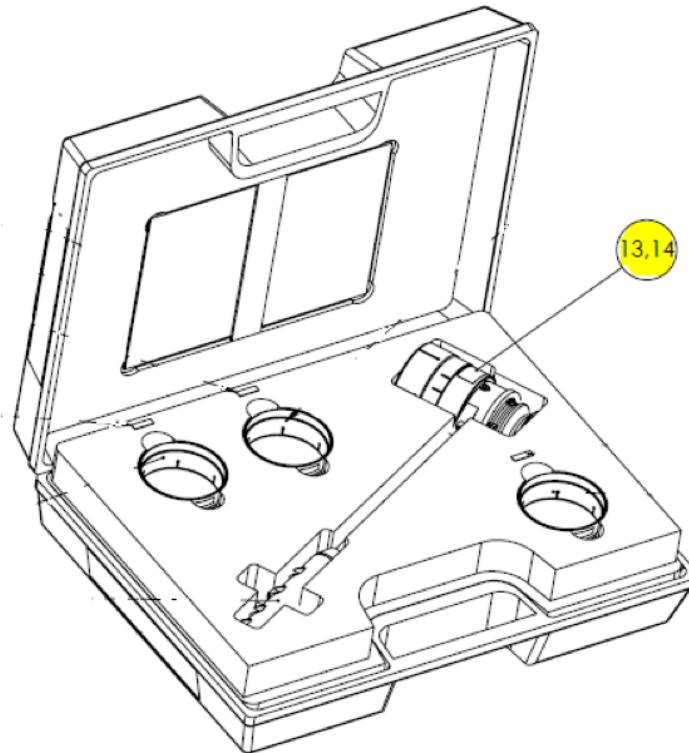

7	RM00730-02	CE mark + European Authorized Representative label	RD006-0205	1
6	RM00843-02	ETS 2.8X-FF Light Sticker on the Transportation Box (EU)	RD005-0302	1
5	US: RM00762-02/ EU: RM00874-02	ETS 3.0X FF Case Identification Label	US: RD006-0202/EU: RD006-0204	1
3	RM00274-02	Fragile	-	1

**SOP**

**Title: Labels Specifications for ETS 2.8X-FF Light (EU)**

Reference code: 100120006

Revision.0

